



**Crozz**  
*two*



# Product specifications

## Crozz series

The Crozz series product line is drawn up by industrial designers in close co-operation with the hospital specialists and suppliers of medical equipment. This, combined with our long time experience in producing sturdy metal construction furniture, provides this new product line the higher level of function, flexibility and design demanded by today's healthcare professionals.

Our products are made from the best quality materials that meet the highest demands for medical and public use. The Crozz series line is based on our popular modular concept, providing you the parts to assemble a cart for any application


## Crozz-two



The Crozz Two is designed for a compact installation of medical (endoscopic) equipment. Based on the sizes of common equipment, this cart is available in three shelf sizes: 320mm, 430mm and 650mm wide.

## Crozz-two basic cart

De Crozz Two is a compact and multi functional equipment cart. Due to a series of modular parts, the Crozz Two can be assembled according to specific wishes and requirements. The basic cart is shown below.

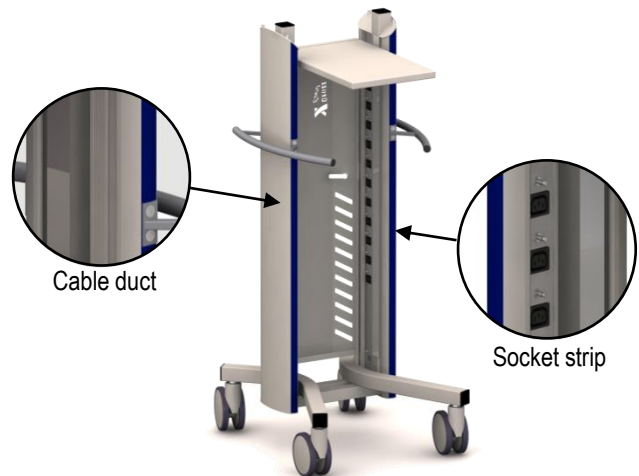

## Characteristics

Sizes wxdxh	60/72/89 x 70 x 142	[cm]
Weight	46	[kg]
Electrical characteristics		
- Rate voltage	230	[VAC]
- Max output power	1600	[VA]
Classification	EN 60601-1	
Earth potential balancing	POAG-6	

## Features

The Crozz Two can be assembled according to your specific wishes and requirements, such as:

- ✓ Shelves, drawers, scope/serum holders etc.
- ✓ Medical isolation transformer
- ✓ Lift and swivel arms for monitor mounting
- ✓ Protection cover for monitors.

**JANSEN**  
MEDICARS

Postbus 49, 3600 AA Maarssen, Nederland  
tel. +31 30 2613500 fax. +31 30 2616741  
url. [www.medicars.com](http://www.medicars.com)